

Shade Matching in Restorative Dentistry: The Art and Science

**Dr Bharat Bhushan¹, Prof. (Dr) Harpreet Singh²,
Dr Gursandeep Kaur Sandhu³, Dr Mohit Chaudhari⁴,
Dr Shubham Sharma⁵, Dr Subreena Zahoor⁶, Dr Sandeep Kumar⁷**

¹Postgraduate Student, Desh Bhagat Dental College and University, Mandi Gobindgarh, Punjab, India

²Dean Research, Professor and Head of the Department, Desh Bhagat Dental College and University, Mandi Gobindgarh, Punjab, India

³Professor, Desh Bhagat Dental College and University, Mandi Gobindgarh, Punjab, India

^{4,5,6,7}Postgraduate Student, Desh Bhagat Dental College and University, Mandi Gobindgarh, Punjab, India

ABSTRACT

Esthetic dentistry seeks to achieve harmony between illusion and reality by creating smiles that integrate teeth, gums, and facial features naturally. A central aspect of this process is color perception, which depends on the interaction of light with dental materials and the human eye's interpretation. Color in dentistry is defined by three primary attributes—hue, value, and chroma—described by the Munsell color system, while additional optical properties such as translucency, opalescence, and fluorescence significantly influence restorative outcomes. To standardize shade matching, both visual and instrumental methods are employed. As esthetics becomes increasingly important to patients, mastering shade selection and understanding color science are essential for predictable, natural, and harmonious dental restorations.

KEYWORDS: Colour, Chroma, Hue, Munsell system, Value

INTRODUCTION

Achieving the right balance between illusion and reality lies at the heart of esthetic dentistry. This harmony is referred to as perception—the aspect related to how appearance is visually interpreted. The primary goal is to create an attractive smile—not just by shaping beautiful teeth, but by ensuring the teeth have natural proportions and blend seamlessly with the gums, lips, and face⁴.

Effective esthetic restorations rely heavily on understanding how light interacts with materials and how the human eye and brain perceive color. Over recent decades, the importance of color research in dentistry has grown significantly. Ultimately, the goal of esthetic dental work is to produce a smile that is both visually appealing and harmoniously integrated. This involves aligning the teeth, gums, and facial features in a natural balance. This article explores the essential color properties that must be considered in esthetic dentistry, along with various tools used for shade selection.

NATURE OF COLOUR:

Color is explained by making reference to a sensation which is captured by our eyes. The human eye is an organ that is receptor of images obtained from an electromagnetic radiation that we refer to as light. It actually corresponds to a narrow segment between the 400 and 800 nm wavelengths approximately of the entire spectrum. We call it as rainbow. When white light is illuminated to object, the color we see has not been absorbed and that light is reflected to our eyes is color is object. So this fact tells us that the quality of light plays an important role in perception of color of an object ¹

MEASUREMENT OF COLOUR

Color is usually described according to the Munsell color space in terms of hue, value, and chroma (Munsell Color System) (Fig.1). When color is determined using the Munsell system, value is determined first followed by chroma. Hue is determined last by matching with shade tabs of the value and chroma already determined.²

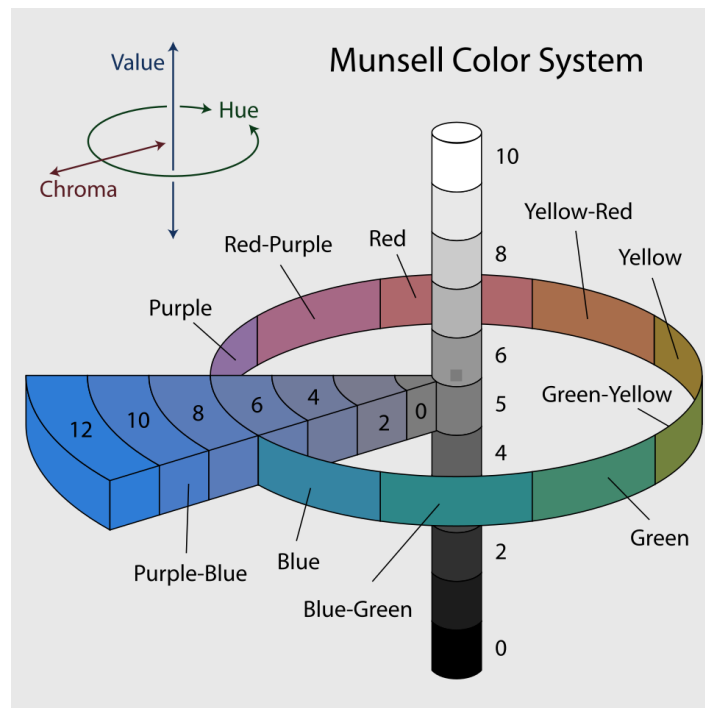


Fig.1 Munsell Color System

Hue

It is the quality that distinguishes one family of color from another. It is specified as the dominant range of wavelengths in the visible spectrum that yields the perceived color, even though the exact wavelength of the perceived color may not be present . Hue is a physiologic and psychologic interpretation of a sum of wavelengths. Hue is represented by A, B, C or D on the commonly used Vita Classic shade guide.

Value

“Value” or brightness, is the amount of light returned from an object . Munsell described value as a white-to-black gray scale. Bright objects have lower amounts of gray and low-value objects have larger amounts of gray and will appear darker. The brightness of a prosthesis is usually increased in two ways: by lowering chroma or by increasing the reflectivity of the surface. Lowering value means less light ret-

urns from the illuminated object and the remaining light is being absorbed or scattered elsewhere.

Chroma

“Chroma” is the saturation, intensity or strength of the Hue. If any dye (say red) is added into a glass of water and the same dye is added again and again, the intensity increases, but the color remains the same (hue). As more dye is added, the mixture appears darker; thus, the increase in chroma has a corresponding change in value. As chroma is increased, the value is decreased; chroma and value are inversely related. Higher numbers on the Vita Classic shade guide represent increased chroma.⁶

In addition to three primary color attributes such as lightness, hue, and chroma, many other optical properties of tooth, such as translucency, opacity, opalescence, surface gloss and fluorescence, should be considered.

FACTORS AFFECTING TOOTH SHADE

Translucency

Translucency is also one of the fundamental factors influencing the esthetic performance of dental restorations. Translucency is the relative amount of light transmission or diffuse reflection from a substrate surface through a turbid medium. Human teeth are characterized by varying degrees of translucency, which varies between transparent and opaque. Increasing the translucency of a crown lowers its value because less light returns to the eye³. With increased translucency, light is able to go by the surface and is scattered within the body of porcelain². The translucency of enamel is influenced by factors such as the angle of incidence, surface texture, polish, wavelength, and dryness level.⁵

Opalescence

Opalescence is the fact in which a material appears to be one color when light is reflected from it and another color when light is transmitted through it. This effect is most commonly observed in enamel and contributes to enhancing the brightness, vibrancy, and depth perception of teeth. Opalescent effect is also called as light scattering effect. It depends on refractive index, size and quantity of opalizers⁶.

Flourescence

Fluorescence is the absorption of light by a material and the spontaneous emission of light in a longer wavelength. UV light is absorbed and fluoresced back as light primarily in the blue end of the spectrum, but it will occur at all wavelengths. The more the dentin fluoresces, the lower the chroma. Fluorescent powders are added to crowns to increase the quantity of light returned back to the viewer, block out discolorations, and decrease chroma².

Surface texture

Young teeth exhibit great surface characterization including stippling, ridges, striation and developmental lobes evidences. The texture of natural teeth is made up of a set of larger and smaller surface fluctuations, which will undergo age related alteration, a considerable impact on reflection and hence on tooth color.⁵

CIE LAB COLOR SYSTEM

The Commission International de l’Eclairage (CIE), in 1931, introduced a system of instrumental color measurement for describing colors. One of the first color spaces defined is CIE XYZ. One of the first clinical implications of the system lies in the development of the VITA 3D Master shade guide

The CIE color system represents a uniform color space, with equal distances corresponding to equal perceived color variances. In this three-dimensional color space, the three-axis are L*, a*, and b*. The

L* represents the value which is a lightness measurement of an object and is quantified on a scale as a 100 L* of value for a perfect black. While the a* value is a quantity of greenness (negative a*) or redness (positive a*). Whereas the b* value is a quantity of blueness (negative b*) or yellowness (positive b*) The a* and b* coordinates approach zero for natural colors and increase in degree for more intense or saturated colors. CIE Lab system’s main advantage is that color variances could be expressed in units that can be related to visual perception and clinical importance¹

COLOR MEASUREMENT

There are two main methods for color matching, the visual method using shade guide, and the instrumental method using color matching devices.

Conventional Visual Method: Visual color matching involves a direct visual comparison of the different color samples in a shade guide with natural teeth in order to determine which tab, or combination of tabs, constitutes the best match .The Munsell color system is a common system for visual color assessment. It is economical, commonly available, and it efficiently compares teeth color with a standardized reference shade guide. The shade guide’s main advantage is that generally they can be used almost anywhere and less expensive. The most popular shade guides currently used are: There are 3 main shade guides used in the daily practice, namely the ones developed by VITA, Bioform by DeTrey and Chromascop by Ivoclar. They are based on Munsell’s principles, which divide the color space into three dimensions – name of the color (hue), the density of the color (Chroma), and vitality (value)⁷

VITA Classic (Fig.2): The color space in Vita Classical is divided into four groups – A, B, C and D based on the dominating HUE (name of the color). In group “A” these are red and brown, in group B – red and yellow, in group C – grey and in group D red and grey. Each of the “letter” groups has subdivisions indexed with Arabic number⁴. With the increase of the number, the Chroma increases, while at the same time, the value decreases.



Fig.2 VITA Clascal Shade Guide

Alternative tab order according to value is B1, A1, B2, D2, A2, C1, C2, D3, A3, D4, B3, A3.5, B4, C3, A4, and C4.

VITA 3D Master (Fig.3) : This shade guide consists of 26 tabs arranged in 5 groups according to their value. Inside the groups, the shade tabs are ordered along two axes – vertical according to Chroma and horizontal based on hue. The first group has two tabs, the second, third and fourth have seven each, while the fifth has three. Each shade tab has three labels – a number ranging from 1 to 5, showing the group and value level, which decreases as the number increases. The letter M defines a mean color tone, which is representative for the given value, while R and L are arranged relative to the median color estimate. The last numeric indicator – 1, 1.5, 2, 2.5 and 3 shows the Chroma levels that are increasing incrementally. Three additional tabs are added for bleached teeth – 0M1, 0M2 and 0M3, the latter denoting the value of 0, three levels of Chroma and mean color tone (M).¹⁴



Fig. 3 VITA 3 D Master Shade Guide

Chromascop (Ivoclar Vivadent, Amherst, NY, USA) (Fig.4): In this shade guide, the color space is distributed in five groups based on the Hue dimension. The color-coding is as: In group 100 – white, In group 200 – yellow, In group 300 – grey, In group 400 – grey, In group 500 – dark brown



Fig.4 Chromascope Shade Guide

Each group has five shade tabs. With the increase in number the chroma of the color increases and at the same time the value decreases. The Shade guide consists of 20 tabs. The number for each of them is a combination of the subgroup and the other two color dimensions. As an addition, accessory shade tabs are added for bleached teeth ¹² – 010, 020, 030 and 040

Shade Guide Made From Composite Material (Fig.5): Most of the manufacturers support the claim that the color of their composite resin or acrylic materials correspond to the ceramic color standards.¹¹

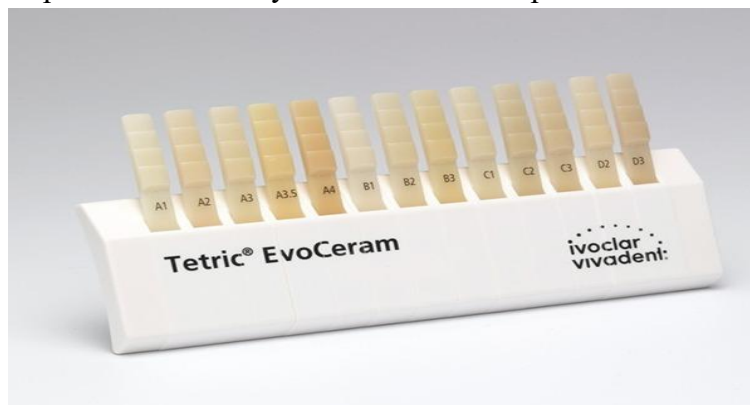


Fig.5 Tetric Evo Ceram Shade Guide

Examples : Esthetic-X (Dentsply), TetricEvo Ceram (Ivoclar) and Venus(Hereus Kulzer)

Instrumental Technique:

In this system, the color space consists of three coordinates: L*, a* and b* . The L* refers to the lightness coordinate, and its value ranges from 0 for perfect black to 100 for perfect white. The a* and b* are the chromaticity coordinates in the red–green axis and yellow–blue axis, respectively. Positive a* values reflect the red color range and negative values indicate the green color range. Similarly, positive b* values indicate the yellow color range while negative values indicate the blue color range.

Shade taking devices: Color-measuring devices usually consist of a detector, signal conditioner and software that process the signal in a manner that makes the data usable in the dental operator or laboratory.

a) Colorimeters

Filter colorimeters generally use three or four silicon photodiodes that have spectral correction filters. These filters act as analog function generators that limit the spectral characteristics of the light striking the detector surface. Shade Eye is an example of a colorimeter based on the natural color concept.²

b) Digital camera as filter colorimeter

The digital camera technology is also being used for shade matching. Instead of focusing light on the film to create a chemical reaction, digital cameras capture images using charged coupled devices (CCDs), which contain millions of microscopically small light-sensitive elements. To get a full color image, most sensors use filtering to look at the light in its three primary colors in a manner analogous to the filtered colorimeter. There are several ways of recording the three colors in a digital camera. The highest-quality cameras use three separate sensors, each with a different filter over it. Light is directed to the different filter/sensor combinations by placing a beam splitter in the camera. The beam splitter allows each detector to see the image simultaneously. The advantage of this method is that the camera records each of the three colors at each pixel location¹¹.

The ShadeRite Dental Vision System and ShadeScan combine digital color analysis with colorimetric analysis, but SpectroShade is the only one that combines digital color imaging with spectrophotometric analysis.

c) Spectrophotometers and Spectroradiometers:

These are instruments designed to produce the most accurate color measurements. Spectrophotometers differ from spectroradiometers primarily because they include a stable light source. There are two types of basic designs commonly used for these instruments. The traditional scanning instrument consists of a single photodiode detector that records the amount of light at each wavelength. The light is divided into small wavelength intervals by passing it through a monochromator. A more recent design uses a diode array with a dedicated element for each wavelength. This design allows for the simultaneous integration of all wavelengths. Both designs are considerably slower than filter colorimeters. However, these are the routinely used color-measuring devices. Eg. VITA Easy Shade (Fig.6), Crystal Eye⁸. Some other available devices are: Shofu Shade Chroma meter, The Shade T Scan, Shade Rite Dental Vision System, The Spectro Shade.



Fig. 6 VITA Easy Shade

CONCLUSION:

Accurate shade matching is essential for achieving natural-looking and esthetic dental restorations. It requires careful evaluation of color under proper lighting, consideration of patient-specific factors, and the use of reliable shade guides or digital tools. Precise shade matching enhances patient satisfaction and the overall success of restorative treatment.

References:

1. Vadher R, Parmar G, Kanodia S, Chaudhary A, Kaur M, Savadhariya T. Basics of color in dentistry: A review. *IOSR J Dent Med Sci*. 2014 Sep;13(9):78-85
2. Johnston WM, Reisbick MH. Color and translucency changes during and after curing of esthetic restorative materials. *Dental Materials*. 1997 Mar 1;13(2):89-97.
3. Todorov R, Peev T, Zlatev S. Shade guides used in the dental practice. *Journal of IMAB–Annual Proceeding Scientific Papers*. 2020 Jun 1;26(2):3168-73..
4. K ro glu A, Makhloota M, Bal BT. A review of color matching in dentistry. *Medical Records*. 2021;3(1):44-9.
5. .Jadhav SD, Pandit VS, Gaikwad AA, Ingale VP, Patil RP, Kaul MS. Shade Selection, a Key to Esthetics. *Journal of the International Clinical Dental Research Organization*. 2024 Jul 1;16(2):110-4.
6. Rajan N, Krishna SR, Rajan A, Singh G, Jindal L. Shade Selection–Basic for Esthetic Dentistry: Literature Review. *Int J Contemp Res Rev*. 2020;11(09):20863-8.
7. .Thejashree N, Vanamala N, Prasad BK. The Importance of Color in Dentistry: A Review. *Asian J. Den. Sci*. 2023;6(1):376-84.
8. Anand D, Kumar GP, Anand DY, Sundar MK, Sharma R, Gaurav A. Shade selection: spectrophotometer vs digital camera–a comparative in-vitro study. *Ann Prosthodont Restor Dent*. 2016 Jul;2(3):73-8.
9. Sproull RC. Color matching in dentistry. Part I. The three-dimensional nature of color. *Journal of Prosthetic Dentistry*. 2001 Nov 1;86(5):453-7.
10. Paolone G, Orsini G, Manauta J, Devoto W, Putignano A. Composite shade guides and color matching. *Int J Esthet Dent*. 2014 Jan 1;9(2):164-82.
11. Tam WK, Lee HJ. Dental shade matching using a digital camera. *Journal of dentistry*. 2012 Dec 1;40:e3-10.